

# The Lojer Capre MC/RC

The Lojer Capre chair range offers a unique, highest quality solution for recovery after day surgery, dialysis and cytostatic treatment, and general examinations. The versatile Capre chairs help bring an air of quality to any facility and improve the patients' comfort during treatment.

The control panel allows adjusting the chair's height, back section, seating angle and both leg sections, independently of one another. The head and leg sections also feature length adjustment. The chair features three preprogrammed positions (incl. beach chair and full trendelenburg position) as well as a memory space for the customers to program and store their own positions.

The Capre chair has an optional multifunction armrest that can be rotated 180 degrees with gas spring assisted angle adjustment and IV-pole attachment. For the foot section there are three options: motorised one-piece section or split leg sections with motors or gas springs. The patient's comfort is further enhanced by an anatomic dual coloured SKAI upholstery, an iPad stand with USB charger, a reading light as well as a coffee tray. If necessary, the chair can also be adjusted into a horizontal position.

*"Highest quality solution for recovery after day surgery, dialysis and cytostatic treatment, and general examinations."*



Capre MC Medical Chair



Capre RC Recovery Chair



Scan the QR code to see the product video, or visit <http://youtu.be/EbrNPcGUhyk>

**LOJER**® For easy care

■ Healthcare | Senior Care | Medical Wellness

## Technical specifications


- **Total length:** 1860...2130 mm
- **Height adjustment:** 600...850 mm
- **Trendelenburg:** 12°
- **Head support adjustment:** 0°...10°
- **Back section adjustment:** -12°...+90°
- **Back section height:** 600 mm
- **Seat section adjustment:** +12°...0°
- **Seat section width:** 620 mm
- **Leg section adjustment:** -90°...0°
- **Armrest adjustment:** 0°...180° rotation, -45°...+20° angle
- **Safe working load:** 160 kg
- **Castors:** Ø 75 mm, central locking
- **USB-output, 4 W**
- **Fulfils MED 42/EEC**
- **Country of origin:** Finland



Features	Capre MC Medical Chair	Capre RC Recovery Chair
Anatomic design upholstery	x	x
Split leg sections	-	x
One-piece leg section	x	-
Basic arm support with rotation	x	-
Versatile arm support with rotation and angle adjustment	-	x
Five electric motors	x	x
Intelligent control system with preprogrammed positions	x	x
Design castors, central locking	x	x
<b>Optional features:</b>		
Foot support for one-piece foot section	o	-
Versatile arm support	o	x
Head support with breathing hole	o	o

Explanation of markings: x = standard feature, o = optional accessory, - = not available

## Accessories

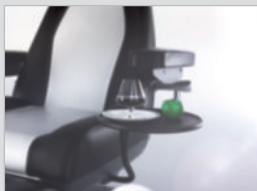

Coffee tray




Tablet support  
with flexible arm



Stitching



Led light

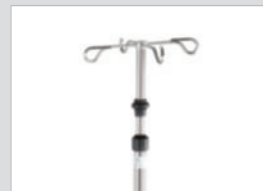

Infusion stand