

Modux Folding Care Bed

Best-selling folding care bed



Ready for use in just a few minutes, without tools



Main features:

1. Small footprint for storage (just 60 cm x 90/100 cm when folded for transport or storage)
2. Folds using electric motors at the press of a button, no tools required
3. Available in 80 cm and 90 cm widths

The patented Modux care bed folds quickly and easily for transport or storage. The Modux is the perfect care bed for units loaning equipment for home care, or as a flexible solution for nursing homes. Although the Modux is a folding bed it still has all the same features as other nursing beds, including a very low minimum height (from 27 cm). Patented Modux nursing beds are flexible and versatile, providing a pleasant care experience for the patient and ergonomic working for care staff. The Modux can also be adjusted to a comfortable position similar to a beach chair.

LOJER® For easy care


Technical data	
Transportation dimensions (w x l)	60 x 90 cm
Overall width with rails	90 or 100 cm
Overall length	224 cm
Height adjustment range	29 - 69 cm (125 mm castors) 27 - 67 cm (100 mm castors)
Recommended mattress size	205 x 78/88 x 13 cm
Total weight	92 or 96 kg
Castors	Ø100 or Ø125 mm
Safe Working Load (SWL):	170 kg
Electrical system	24V
Back section adjustment	70°
Foot section adjustment	38°
Middle section adjustment	40°
Anti-trendelenburg	12°
CE -marked	Yes
Warranty	2/10 years
Country of manufacture	Finland

Factory installed options:

- Ø100 mm castors

Other accessories:

- Lifting pole with triangle handle
- Lift support
- IV-pole with one hand quick adjustment
- Reading and dining tray, on top of the side rails
- Reading lamp
- Padded coating on top of the side rail, one unit
- Bed ends: birch

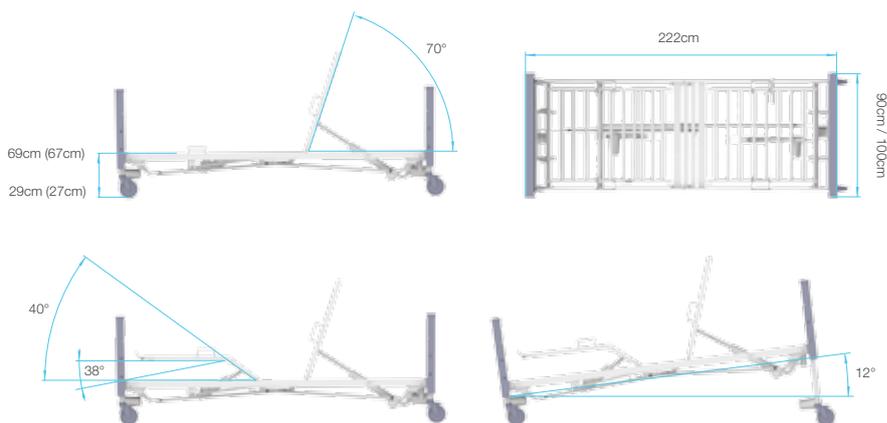

LOJER® For easy care