

ScanAfia XHS Nursing Bed



Product pictured with optional accessories

Main features:

1. High Safe Working Load (SWL) of 260 kg
2. Low minimum height (from 32 cm) improves patient safety and self-reliance
3. Back section sliding mechanism prevents patient slipping downwards and improves comfort

The ScanAfia XHS nursing bed is the number one quality choice for the healthcare professional. The newest nursing bed in the product range, it has been designed to suit a variety of care environments, where design and appearance are appreciated as much as reliability and ease of use. Although it is able to take up to 260 kg weight (SWL), the ScanAfia XHS is lightweight and easy to manoeuvre, thanks to the 125 mm castors fitted as standard.

The bed's exceptional height adjustment range makes it safer and more comfortable to use for recuperating patients and makes work easier for carers. Our many years' experience in designing and manufacturing hospital and nursing beds ensures that the ScanAfia XHS is made for long-term use with minimal maintenance. The base of the bed is made of hygienic ABS-plastic with ventilation holes. This base structure allows the foam mattress to retain its characteristics better than it would with a traditional mesh base. The base can easily be removed for disinfecting.

The Range:

- ScanAfia XHS 280/290: 2-part bed section (205 X 80/90 cm), with electrical height and back section adjustment
- ScanAfia XHS 480/490: 4-part bed section (205 X 80/90 cm), with electrical height, back, calf and leg section adjustment

LOJER® *For easy care*


Technical data and features	
Total width (incl. siderails)	92 cm / 102 cm
Patient surface width	80 cm / 90 cm
Total length	217 cm Long model: 227 cm
Total weight	130 - 132 kg
Height adjustment range	Min. 32 cm - max. 82 cm (depending on the castor size and type)
Recommended mattress size	78-80 cm or 88-90 cm x 205 cm x 9-13 cm
Safe Working Load (SWL)	260 kg
Castors	Central locking Ø125 mm with directional wheel
Back section adjustment	70°
Thigh section adjustment	30° (4-part bed)
Calf section adjustment	11° (4-part bed)
Anti-Trendelenburg (optional)	20°
Other standard features	Hand control bed ends: Classic beech or Classic birch two double-locking wooden side rails with locking button quick release (CPR) for back section back section sliding mechanism removable base with mattress attachments fixing for 2 lifting poles and IV stand fixing for restraining straps vertical raising motion
Electrical system	24V
CE-marked	Yes
Country of Manufacture	Finland

Factory installed options:

- Ø125 / Ø150 mm twin castors
- Anti-Trendelenburg
- Battery driven (short-term backup)
- Set of bumper rollers
- Scandic white oak bed end
- Special length (extendable) model

Other accessories:

- Lifting pole with triangle handle
- Lift support
- IV pole with one hand quick adjustment
- Reading lamp
- Padded side rail (one unit)
- Urine bag holder
- Reading and dining tray on top of the side rails
- Nurse locking panel (allows functions to be locked to prevent patient use)
- Flexible hand control holder

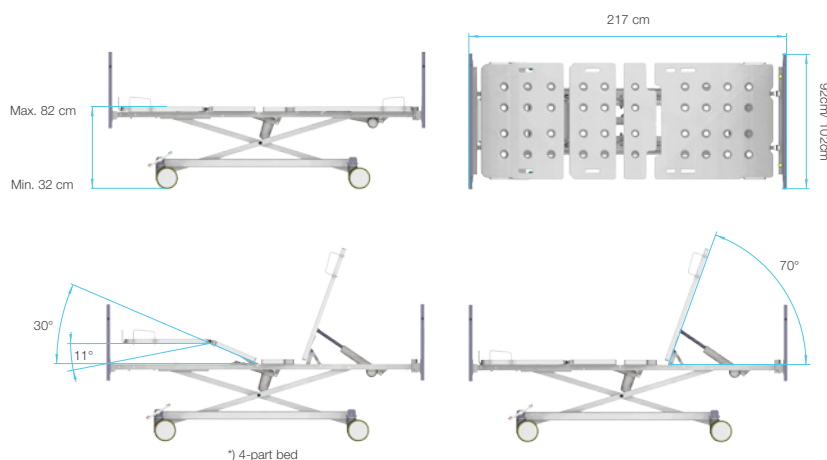

LOJER® For easy care