



MERIVAARA KIDS MATTRESS

Hygienic mattress

Baby and Junior mattresses

Merivaara codes 902 / 905 / 908

Purpose of use


- Good basic mattresses for newborns and kids.

Product structure

- Core 902 & 905 mattresses R25085 (25 kg/m³) fire retardant HR-foam.
- Core 908 mattress D33145 (33 kg/m³) fire retardant foam.
- Cover polyurethane coated hygienic bi-elastic fabric, which is impermeable to liquids but still breathable. Material is bi-elastic and treated against bacteria and fungi (active biocidal substance is zinc pyrithione (EU 528/2012). Zipper on the end.

Benefits

- Secure and easy-care mattresses can be used on both sides.
- Facilitate and accelerate the work of staff when cover is easy to clean by wiping.

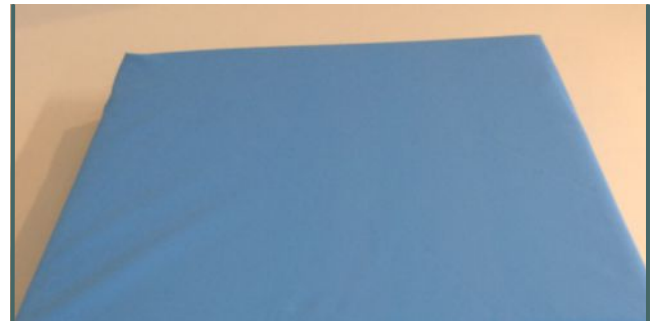

Cleaning and maintenance

- Disinfection by wiping according to the instructions with diluted chlorine-amine or an alcohol-based detergent or ethanol /chlorine max 1000 ppm, rinse and dry thoroughly. If necessary, wash the cover 70-95°C inside out, gentle tumble drying at low temperature.
- Turn the mattress regularly to ensure even wear. Replace broken cover to a new one.
- The bed frame should be lathed or perforated to allow for proper ventilation.

Flame retardancy & CE

Flammability by standard IMO FTPC part 9. Merivaara Kids mattresses are CE-marked. Country of origin Finland.

Size

Standard sizes 38*62*5, 60*120*8 and 75*160*12 cm

