

Bobath XB Mat Table

Sturdy treatment table for neurological therapy



Main features:

1. Very high safe working load (SWL) 400 kg
2. Patented structure that makes the table sturdy and well suited for heavy use
3. Patient lifts fit easily underneath even when the table is at its lowest height (40 cm)


The renewed Lojer Bobath Mat Treatment Table is designed primarily for neurological therapy and Bobath therapy. The table has a patented structure that makes it very sturdy and well suited even for constant heavy use. It is designed in such a way that the therapist can be on the table at the same time as the patient who is being treated – safe working load (SWL) is up to astonishing 400 kg.

The tabletop is upholstered to a high standard and the sides are strengthened with plastic edging. The upholstery material is extremely durable and antimicrobial fabric that is also easy to clean, waterproof and has excellent disinfectant features. Tabletop's foam thickness is 30mm, which makes the patient surface less elastic and more suitable for examinations and treatments of neurological disorders. Bobath's frame structure allows patient lifts to fit easily underneath it (ground clearance is 15 cm) even when the table is at its lowest height. Extremely low minimum height (40 cm) ensures that even patients with wheelchair are able to get onto the table with ease.

There are two models available: 4010XB has one-part tabletop and 4020XB comes with two-part tabletop and adjustable back section. Both models can be equipped with optional face hole to enhance patient comfort. Both models' motors lower using gravity and in emergency situations - if there's something caught between the mechanism – table rises.

LOJER® For easy care

Technical data and features		
Model	4010XB	4020XB
Sections	one-section	two-section
Total width	100 cm / 120 cm	100 cm / 120 cm
Total length	200 cm	200 cm
Total length with face hole (optional)	207 cm	207 cm
Total weight	145 kg	145 kg
Height adjustment range (electric)	40 - 95 cm	40 - 95 cm
Back section length	-	70 cm
Back section adjustment range	-	0-73°
Trendelenburg / anti-Trendelenburg (optional)	20° / 20°	20° / 20°
Castors	Ø 150 mm central locking twin	Ø 150 mm central locking twin
Other standard features	Plastic edging to protect upholstery hand control with safety lock functions (pressing two buttons for 3 sec) safety feature: auto-lifting in emergency situations	Plastic edging to protect upholstery hand control with safety lock functions (pressing two buttons for 3 sec) safety feature: auto-lifting in emergency situations
Upholstery	Several material and colour options, more information on request	Several material and colour options, more information on request
Safe working load (SWL)	400 kg	400 kg
Ground clearance	15 cm	15 cm
Electric system	24V	24V
CE-marking	Yes	Yes
Warranty	2 years	2 years
Country of Manufacture	Finland	Finland

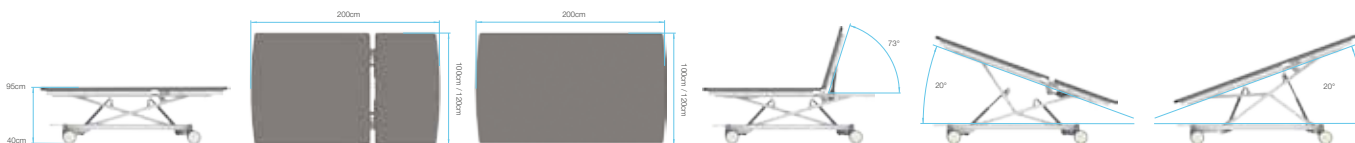

Options

- Trendelenburg and anti-Trendelenburg (incl. battery)
- Foot control
- Face hole
- Battery powered

Accessories

- Accessory rail
- Lift support
- Push handles
- Examination lamp
- Side rails

