

EMERGO™


Lojer Merivaara multipurpose patient trolley



LOJER®

For Easy
Care



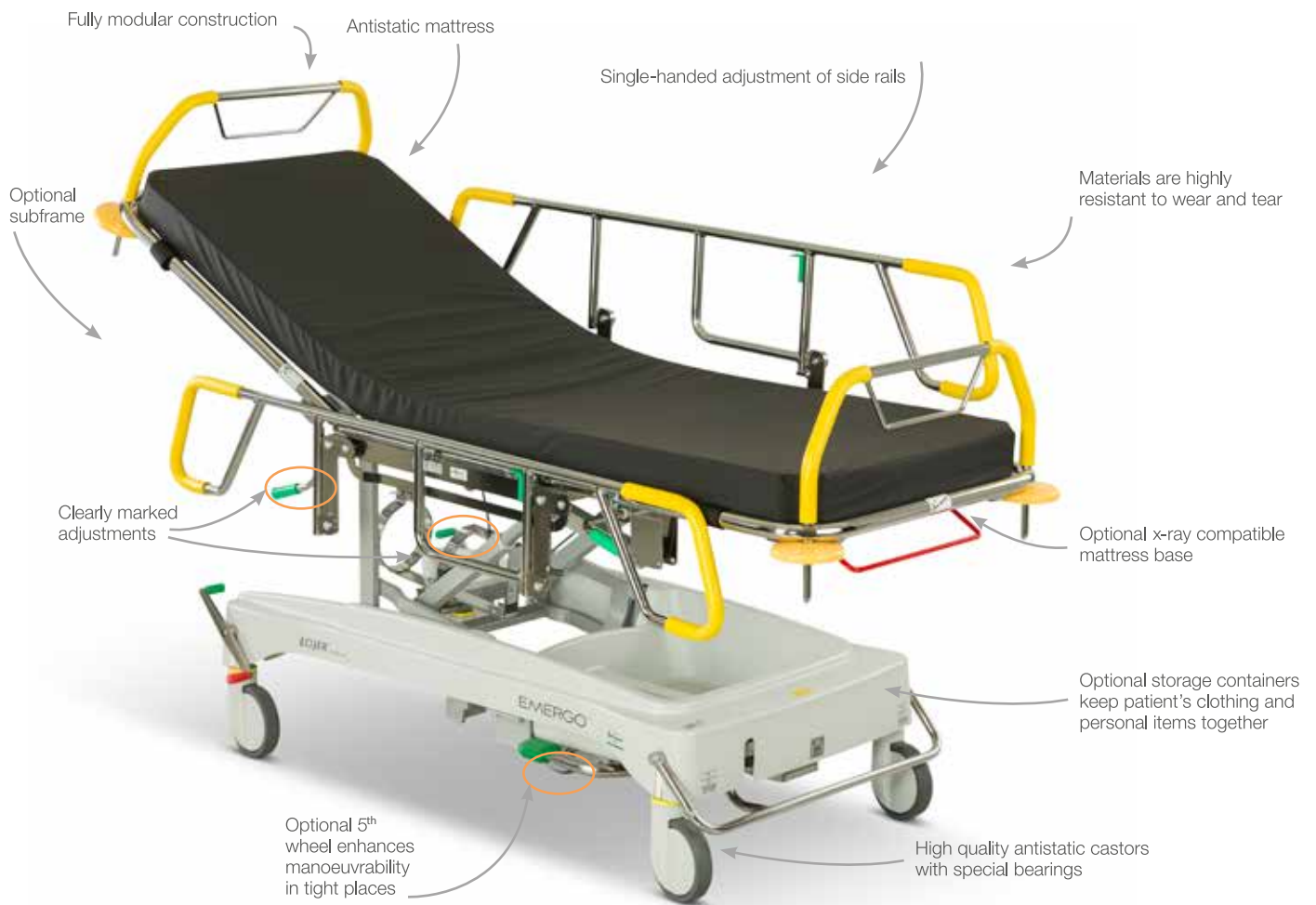

1. Modular structure that allows several different configurations
2. Clear and easily controllable functions
3. It is easy and light to carry out different care positions

The multipurpose Emergo patient trolleys (stretchers) are used for patient transportation, emergency, recovery, ICU and minor procedures. The safe and sturdy construction makes it ideal for heavy use. The versatile range of options and adjustments ensure efficient patient transfer, patient comfort and ergonomic working positions for health-care staff. The adjustable mattress base provides comfortable positions for patients to recover eliminating unnecessary patient transfer from the trolley (e.g. in sitting position).

Large, easy rolling castors provide excellent mobility and the optional fifth wheel improves overall mobility reducing the strain for the staff. Emergo has many attachment points for the IV-rod and other accessories. The wide range of accessories available will serve the demanding needs of nursing staff for many years.

The range:

- 6250: 2-part bed section with hydraulic height adjustment and gas spring assisted Trendelenburg and back section
- 6270: 2-part X-ray bed section with hydraulic height adjustment and gas spring assisted Trendelenburg and back section
- 6280: 2-part X-ray bed section with subframe, hydraulic height adjustment and gas spring assisted Trendelenburg and back section
- 6350: 3-part bed section with hydraulic height adjustment and gas spring assisted Trendelenburg, back section and leg section
- 6370: 3-part X-ray bed section with hydraulic height adjustment and gas spring assisted Trendelenburg, back section and leg section
- 6380: 3-part X-ray bed section with subframe, hydraulic height adjustment and gas spring assisted Trendelenburg, back section and leg section

Factory-assembled options and accessories




Castors central locking antistatic
Ø 200 mm
Fixed, without castors (eg. for ships).



X-ray compatible mattress base,
easy to clean. Available in models
6270, 6280 and 6370



5th wheel for improved mobility

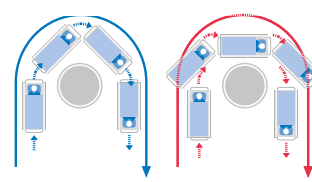

On the left is shown how the turning
radius decreases with the 5th wheel.
On the right the same situation
without the 5th wheel.

Storage



ABS storage box
100006620 large
100006622 large, for dual pedals
100006625 large, for fifth wheel
100006626 small, only with large storage box
(Not compatible with fixed height frame.)



Wire basket under back section
100000582



Wire basket under foot section
100000583
(Not compatible with 5th wheel.)



Accessory basket (requires trolley end)
128009784

Accessory holders and rails



Paper roll holder
100000733



Accessory rail (pair) 700 (replaces side
rail, arm rest in picture not included)
100000700
Arm rest 100000125



Oxygen bottle holder 589
100000589 for 5 kg bottles
100005860 for 2 kg bottles



Infusion rod holder 5785
100005785

Infusion rods




Infusion rod, 4-hooks
100000499



Fixed, foldable infusion rod for
subframe
100006678

Other accessories

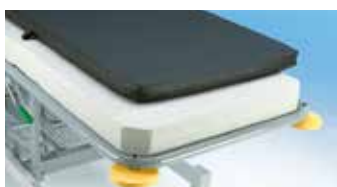

Hygienic and antistatic mattresses
Size 650mm x 1950mm
Thickness 50mm, 80mm or visco-elastic
100mm




Patient restraining straps (set)
100000790



Side rails (pair)
100006713 chromed/yellow



Push bar, foldable
100006793 chromed/yellow
padding
100006795 epoxy coated/yellow
padding

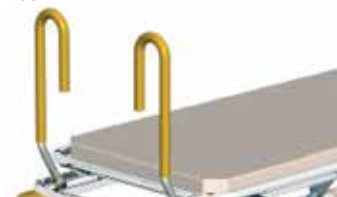

Long push bars, detachable (pair) 6770
10000 6770 chromed, yellow padding



Trolley End 567
100000567

Colours in pictures may differ from those currently available.

| Technical data and features | | | | | | |
|-----------------------------|--|------------------------------------|------------------------------------|--|--|--|
| Model | 6250 | 6270 | 6280 | 6350 | 6370 | 6380 |
| Intended use | Emergency, transport, minor procedures | Emergency, ICU, minor procedures | Emergency, ICU, minor procedures | Emergency, transport, minor procedures, recovery | Emergency ICU, transport, minor procedures, recovery | Emergency ICU, transport, minor procedures, recovery |
| Base sections | 2 | 2 | 2 | 3 | 3 | 3 |
| X-ray platform | No | Yes | Yes | No | Yes | Yes |
| Subframe for equipment | No | No | Yes | No | No | Yes |
| Width | 750 mm | 750 mm | 750 mm | 750 mm | 750 mm | 750 mm |
| Patient surface width | 705 mm | 705 mm | 705 mm | 705 mm | 705 mm | 705 mm |
| Length | 2075 mm | 2075 mm | 2150 mm | 2075 mm | 2075 mm | 2150 mm |
| Recommended mattress size | 650 x 1950 mm | 650 x 1950 mm | 650 x 1950 mm | 650 x 1950 mm | 650 x 1950 mm | 650 x 1950 mm |
| Height adjustment | 560 - 930 mm | 625 - 995 mm | 625 - 995 mm | 560 - 930 mm | 625 - 995 mm | 625 - 995 mm |
| Safe Working Load (SWL) | 225 kg | 225 kg | 225 kg | 225 kg | 225 kg | 225 kg |
| Trendelenburg | 14° | 14° | 14° | 14° | 14° | 14° |
| Anti-Trendelenburg | 7° | 7° | 7° | 7° | 7° | 7° |
| Back section adjustment | 0 - 70° | 0 - 70° | 0 - 70° | 0 - 70° | 0 - 70° | 0 - 70° |
| Leg section adjustment | - | - | - | 0 - 45° | 0 - 45° | 0 - 45° |
| Castors | Central locking, antistatic 150 mm | Central locking, antistatic 150 mm | Central locking, antistatic 150 mm | Central locking, antistatic 150 mm | Central locking, antistatic 150 mm | Central locking, antistatic 150 mm |
| Country of manufacture | Finland | Finland | Finland | Finland | Finland | Finland |
| CE marked | Yes | Yes | Yes | Yes | Yes | Yes |

Other features:

- Yellow bumpers in each corner
- Pedals for hydraulic height adjustment on both sides
- Trendelenburg and anti-Trendelenburg adjustments
- Antistatic central locking 150 mm castors
- Gas-spring assisted back section
- 2- or 3-part base

