

Screen-Systems

Healthcare | 2025



NEW!

Now with exclusive
illustrations
by Lisa Tegtmeier



ropimex®
since 1976

Ideal for all departments



ropimex® Sound insulation: Acoustic solutions
For clinics and office. Request catalog now!

Telescopic screens,
Telescopic arms RTS/RTI..... 4 - 13

Wall brackets /
Brackets for equipment rails
WH..., GSH... 14 - 15

Mobile screens FG95 16 - 17

Vario rail system RVS 18 - 21

Curtains
RTX, VHC, CQ, Z 22 - 25

ropiCARE Folding curtain..... 26 - 27

Folding screens RFW 28 - 35

Butterfly Folding wall RBF 36 - 39

Butterfly Paravent RBP 40 - 41

Paravents RLP..... 42 - 45

Telescopic cubicle RKLE...,
Shower splash guard RD... 46 - 47

About us / Certificates 48 - 51



Certified according to
DIN EN ISO 9001:2015
DIN EN ISO 14001:2015



ropimex products meet the requirements of the
regulatory specifications for medical devices in accordance with the annex to Declaration of Conformity
EU 2017/745, ropimex® visual protection systems



Please visit
our
virtual showroom



Scan



Telescopic arms

You can't trip over things that hang on the wall. Which is why we created ropimex® telescopic screens, the screen that's only there when it's needed.

After use, simply push the screen together and turn it towards the wall. Your room is back to looking its familiar self.

To a max. 210 cm
Model RTS:
swing-out



To a max. 220 cm
Model RTI:
Fold-down and
swing-out



From 220 cm
Model RTH:
with hooks and
ceiling mounts




Telescopic arms

Optimized in handling, cleaning and disinfection




Smooth telescopic tubes and plastics for one hygiene-friendly surface.

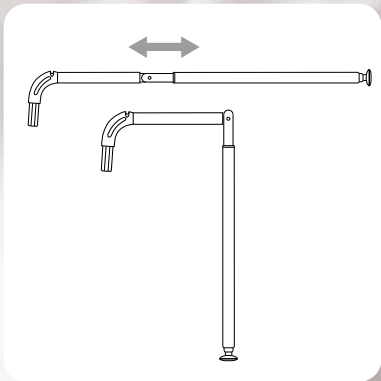

Model RTI: Bending mechanism seamlessly integrated.



If the curtain rings are loaded they release themselves from the eyelets of the curtains and therefore avoid injuries and damage to the system.

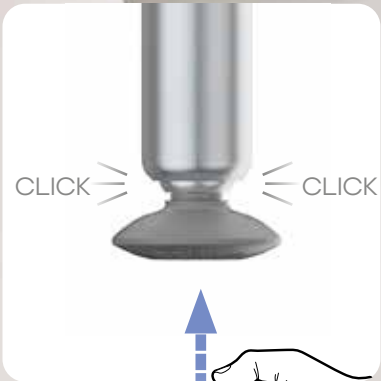

Tail with locking - fixed telescope at rest



Telescopic arms

Model **RTS** - swing-out



Telescopic screen Model RTS - swing-out



ropimex® - telescopic screens are made of anodised aluminium, high-quality, glass-fibre reinforced plastics. Easy to clean and to disinfect.



Model	Telescopic screen extension	Weight	Curtain rings included in the scope of delivery 
RTS 1.2	55 cm → 120 cm 	0,5 kg	
RTS 1.5	65 cm → 150 cm 	0,6 kg	
RTS 2.1	85 cm → 210 cm 	0,7 kg	




IF they are loaded they release themselves from the eyelets of the curtains and therefore avoid injuries and damage to the system.

Wall and rail brackets
see page 14




Recommended curtains see from page 22


RTX ropitex

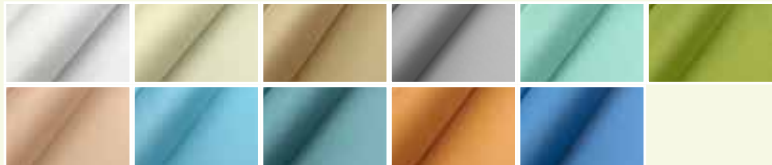

CQ Cleanique



VHC staph CHEK

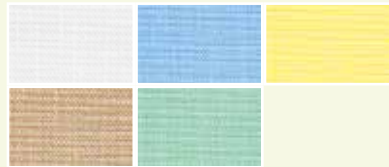

Telescopic arms

Model RTI - Fold-down and swing-out



Curtains

washable

surfaces

disinfectable

with steam

disinfectable


with vaporisation

disinfectable

sustainable

reparable

Telescopic screen Model RTI - Fold-down and swing-out

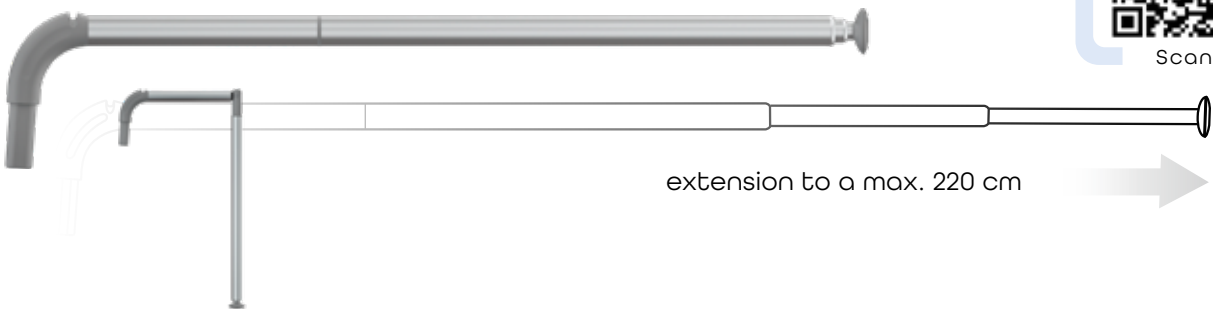

ropimex® - telescopic screens are made of anodised aluminium, high-quality, glass-fibre reinForced plastics. Easy to clean and to disinfect.



Model	Telescopic screen extension	Weight	Curtain rings included in the scope of delivery
RTI 1.3	70 cm → 130 cm	0,6 kg	
RTI 1.6	80 cm → 160 cm	0,7 kg	
RTI 2.2	100 cm → 220 cm	0,8 kg	




If they are loaded they release themselves from the eyelets of the curtains and therefore avoid injuries and damage to the system.

Wall and rail brackets
see page 14




Recommended curtains see from page 22

RTX ropitex

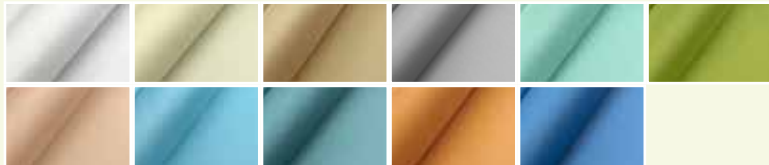

CQ Cleanique




VHC staph CHEK



Telescopic arms

Model **RTH** - Fold-down and swing-out with end support


Used in combination with a height-adjustable ceiling mount RST or with an end support (GL), with the longest telescopic arm RTH 350 you can achieve a screen length of 350 cm.

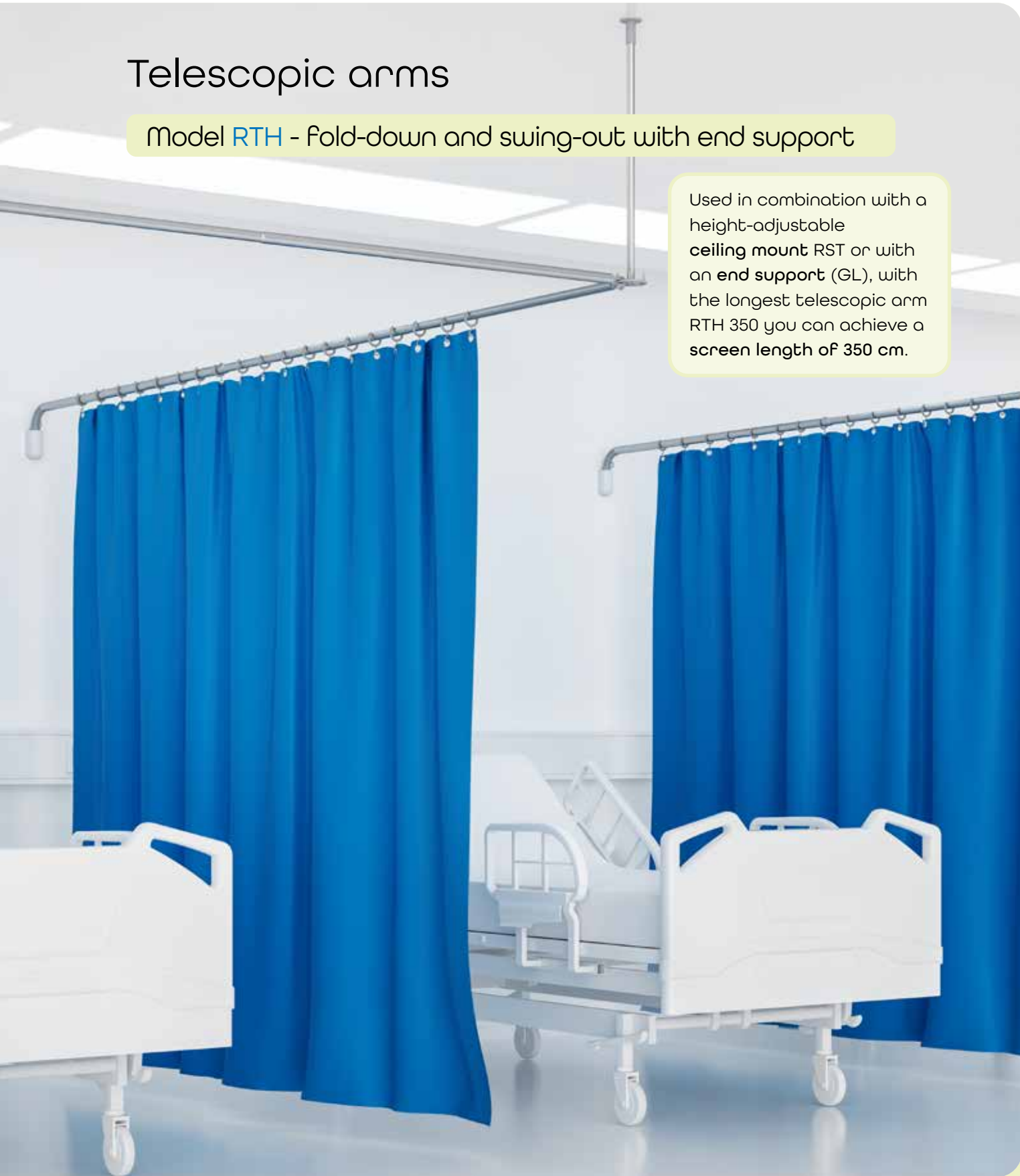

Telescopic screen Model RTH - Fold-down and swing-out with end support



Scan

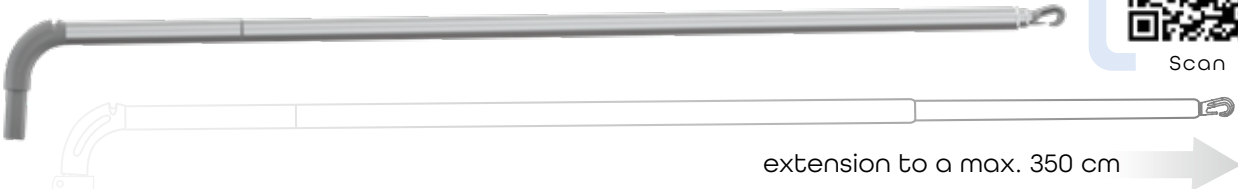

For installation anywhere in the room: ceiling mount



For wall-to-wall installation: end support (GL)



ropimex® - telescopic screens are made of anodised aluminium, high-quality, glass-fibre reinforced plastics. Easy to clean and to disinfect.

Model	Telescopic screen extension	Weight	Included in the scope of delivery
RTH 300	300 cm with hook at the distal end 142 cm → 300 cm	1,0 kg	Wall bracket WH 95
RTH 350	350 cm with hook at the distal end 152 cm → 350 cm	1,1 kg	Ring VR End support GL



Model **RST** - adjustable ceiling mounts

Model	Extension	Room height From - to
RST 4	45 -75 cm	245 - 275 cm
RST 8	75 -135 cm	275 - 335 cm
RST 10*	110 -200 cm	335 - 400 cm

*not illustrated



The glass-fibre reinforced polyamide disc at the bottom end of the ceiling mount is used to hold the joints of the telescopic screens RTH.

Brackets for telescopic screens

Model WH - Wall brackets
Model GSH - Brackets for equipment rails



WH: Wall brackets for telescopic screens



Model		Dimensions
WH 95	Wall bracket for all telescopic screens, die-cast aluminium, concealed screws, adjustable to compensate for unevenness of the wall, colour coated, colours: white (RAL 9010), signal grey (RAL 7004)	Height 10 cm Width 6 cm
WH 95 PL	For walls with low strength, e.g. plasterboard. Wall bracket for all telescopic screens. White or grey, concealed screws, with backplate 20/15 cm	Height 15 cm Width 20 cm
WH	Wall bracket for all telescopic screen	
WH/LA	Wall bracket for bridging service ducts up to 15 cm	Height 16,5 cm Width 6 cm
WH/F	Wall bracket with knurled thumb screw	

GSH: Brackets for equipment rails




Model		Dimensions
GSH/95	Bracket for equipment rails, adjustable tilt	Height 10 cm Width 6 cm
GSH/N	Bracket for equipment rails, wall-supported, height-adjustable	Height 55 cm Width 5 cm
GSH/DK	Bracket with "double clamp" for equipment rails, height-adjustable	Height 55 cm Width 5 cm
GSH/F	Bracket for equipment rail 25/10 mm, fixing by means of clamping screw	Height 13,5 cm Width 4 cm

Mobile screens Model FG95



FG95 trolley

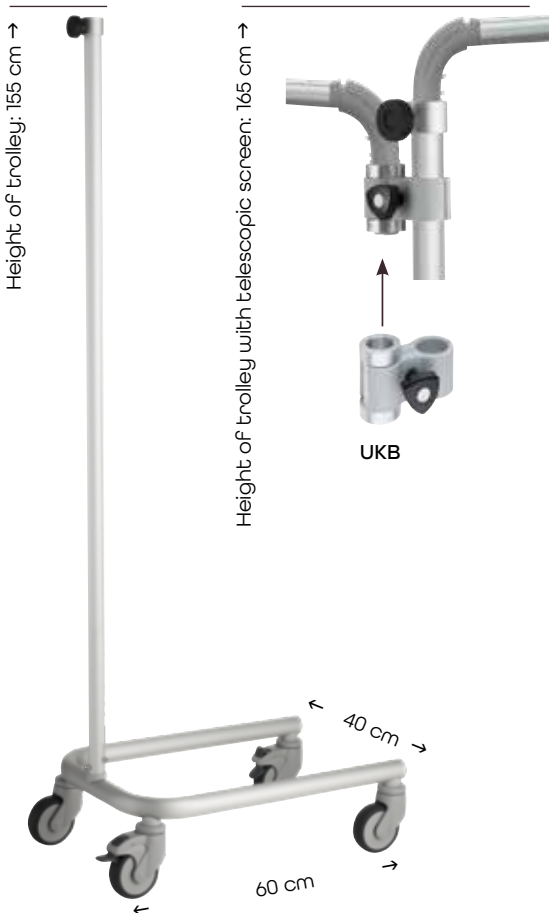

Model

FG95	Trolley , 155x60x40 cm (H/D/W), 4 castors ø 100 mm, 2 of which are lockable. Aluminium anodised, glass-fibre reinforced polyamide, corrosion-resistant, disinfectable. The FG 95 trolley can be fitted with all telescopes of the 'RTI' model series.
LA	Electrically conductive castors available for an extra charge.
UKB	Universal clamp for holding a second telescopic screen - RTI 1.3 or RTI 1.6.

Telescopic screen Model RTI



We recommend the following fitout

With one telescopic screen


RTI 2.2	Telescopic screen, Fold-down, extension 220 cm	CE
---------	--	----

With two telescopic screens and UKB clamp


RTI 2.2	Telescopic screen, Fold-down, extension 220 cm (not suitable for the UKB clamp)	CE
RTI 1.6	Telescopic screen, Fold-down, extension 160 cm	CE
RTI 1.3	Telescopic screen, Fold-down, extension 130 cm	CE

Recommended curtains
see From page 22


RTX ropitex

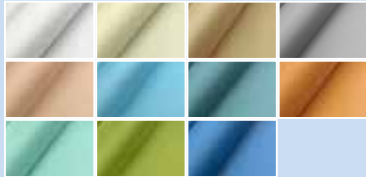


CQ Cleanique



VHC staph CHEK

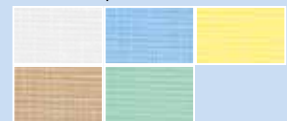

Vario rail system model **RVS**

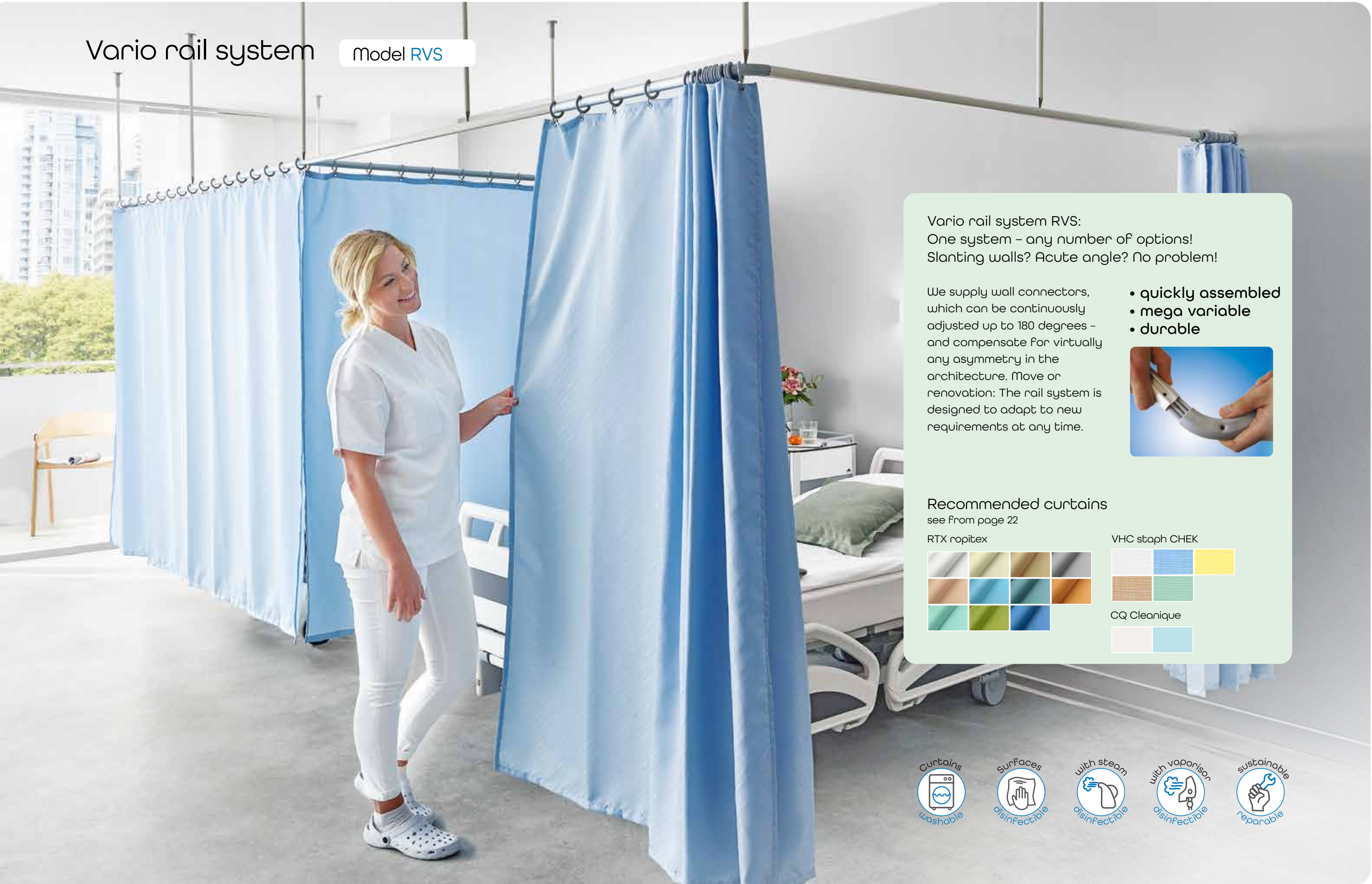

Vario rail system RVS:
One system – any number of options!
Slanting walls? Acute angle? No problem!

We supply wall connectors, which can be continuously adjusted up to 180 degrees – and compensate for virtually any asymmetry in the architecture. Move or renovation: The rail system is designed to adapt to new requirements at any time.

- quickly assembled
- mega variable
- durable




Recommended curtains
see From page 22

RTX ropitex



VHC staph CHEK

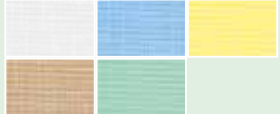

CQ Cleanique



Vario rail system Model **RVS**


Easy ordering: A simple sketch suffices!

Simply prepare a hand-drawn, dimensioned sketch of the required room layout, enter the widths and lengths as well as the room height and the required headroom.

We will put together the complete system for you.



Scan

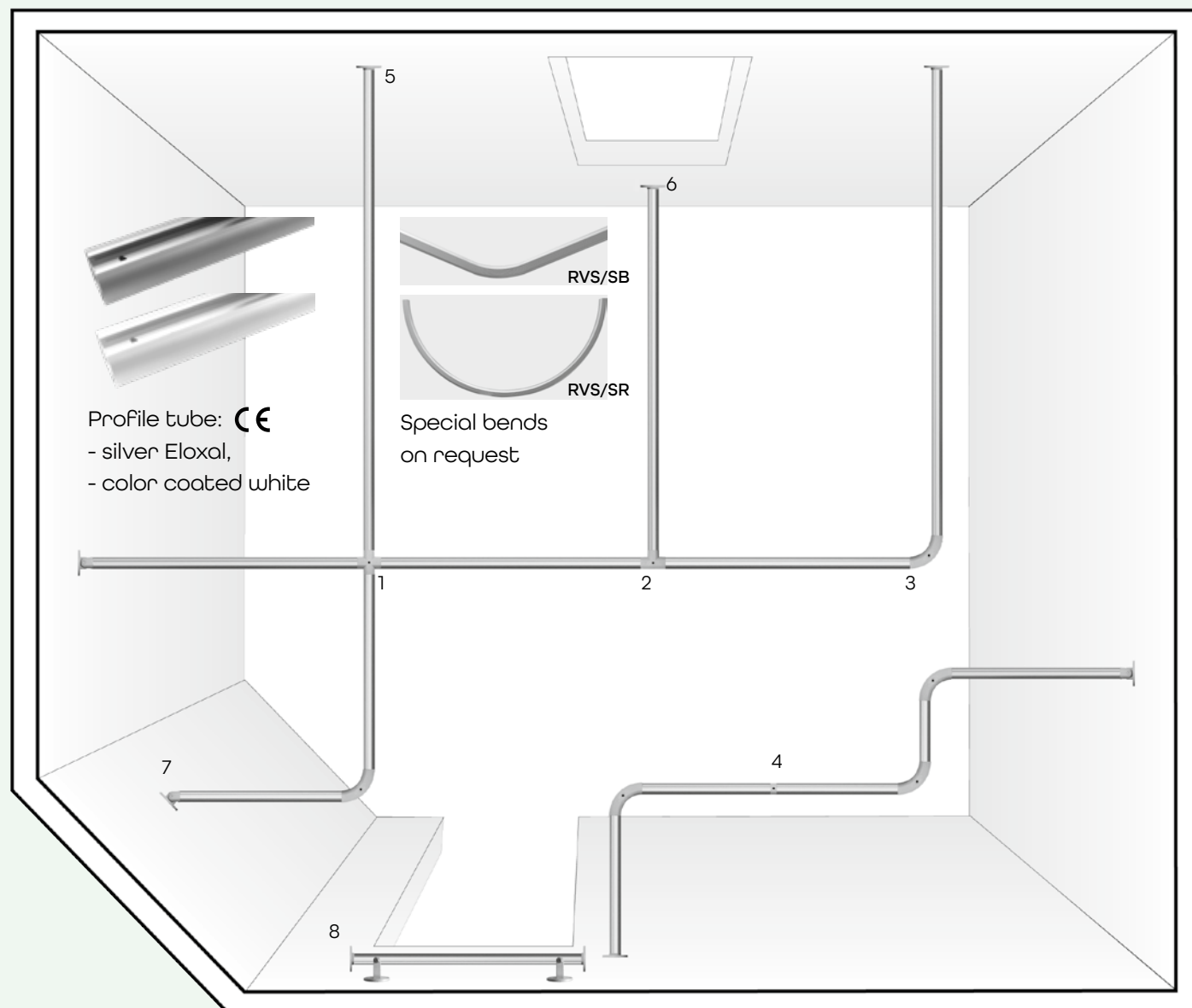

Profile tube: **CE**
- silver Eloxal,
- color coated white

Special bends
on request

RVS/SB


RVS/SR



Curtain ring RVS/VR
made of unbreakable
polyamide: Releases itself
if overloaded (approx.
2 kg).




Ringstopper: For fixation of a curtain ring at the
beginning or the end of the rail.

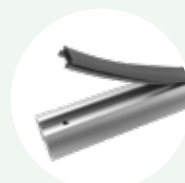


Silicone cover profile
RVS/AP

With the RVS connectors, you can create cabs as you wish

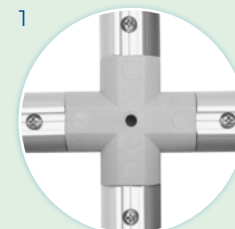


1
Cross-connector
(RVS/KV)

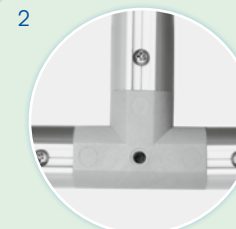

2
T-connector
(RVS/TV)



3
Corner solutions
RVS/BV
(radius 60 mm)

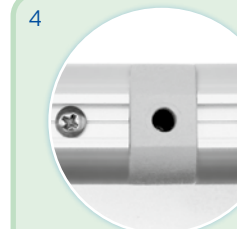

4
Change later:
simply disconnect the
profile tube and ins-
tall the connector of
your choice, e.g. linear
connector RVS/LV


RVS-Wall connections - for every requirement



5
Free-standing end:
End cap RVS/EK
stops the curtain.



6
Straight at the wall:
Wall connector RVS/WA

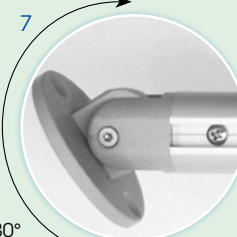

7
Sloped walls, acute angles: the Vario wall connector
RVS/WAV compensates for and levels out asymmetries
up to 180° and length differences up to +/- 5 mm.


+/- 5 mm

180°

Flexible to the ceiling - RVS ceiling connectors



8
**Mounted directly
under the ceiling:**
with ceiling bracket
RVS/DH (distance
from top of ceiling
to bottom of rail
approx. 6 cm).

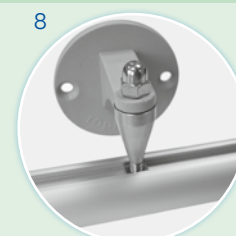

**Easy above the lin-
tel:** e.g. at windows
or doors using the
wall connector
RVS/WAS



**Suspended at
headroom height:**
with the adjustable
ceiling mount
(+/- 20 mm) RVS/DS




**For free-standing
end RVS rails:**
adjustable double
ceiling mount
(+/- 2cm) RVS/DDS



Cover cap for
suspended
ceilings




RVS/DSPL
Reinforcing
disc for ceiling
support for
e.g. plaster-
board ceilings

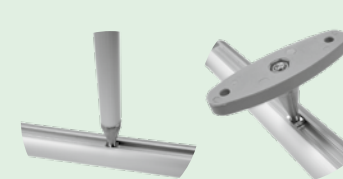

Anchoring on profile rail:
possible along the whole length by
means of square nut; For adjustable
ceiling mount RVS/DS or with ceiling
bracket RVS/DH. Max. spacing: 2 m.



Anchoring via all connectors:
Standard threaded hole for screwing
in the adjustable ceiling mount RVS/DS
or the ceiling bracket RVS/DH
ensures firm support.

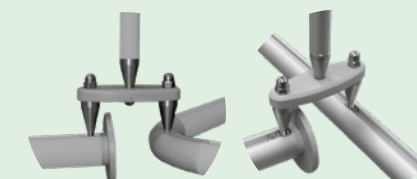

Spacer:
RVS/AH for
Corner- and
T-Solutions.

Curtains



RTX: ropitex BIOACTIVE the new ropimex® curtain for sanitary blocks made of 100 % bioactive polyester fibre



Scan

- Width: 210 cm
- Weight: 140 g/m²



- Care instructions




Standard dimensions

	Width	Height in cm
RTX 110	110	175
RTX 121	140	120
RTX 122	140	140
RTX 123	140	175
RTX 151	170	120
RTX 152	170	140
RTX 153	170	175
RTX 200	200	175
RTX 211	230	120
RTX 212	230	140
RTX 213	230	175
RTX 215	250	175
RTX 217	270	175
RTX 300	300	175
RTX 330	330	175
RTX 360	360	175
RTX 390	390	175
RTX 420	420	175
RTX 460	460	175
RTX 500	500	175

Special dimensions on request.

Available in 11 colours

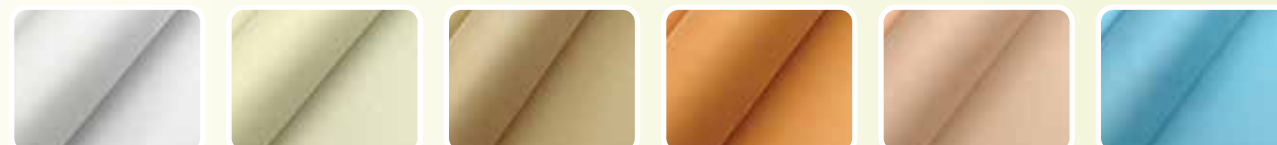

White Champagne Linen Kumquat Peach Light blue



Turquoise blue Ocean blue Pebble Pastel green Pistachio

Flame retardant curtains can save lives!

staph CHEK®

VH/C: ropimex® curtains made of staph CHEK, the fabric for sanitary blocks, showers and bathrooms



Scan

- Width: 180 cm
- Weight: 330 g/m²




- Care instructions



Standard dimensions

	Width	Height in cm
VH/C 121	140	120
VH/C 122	140	140
VH/C 151	170	120
VH/C 152	170	140
VH/C 211	220	120
VH/C 212	220	140

Special dimensions on request.

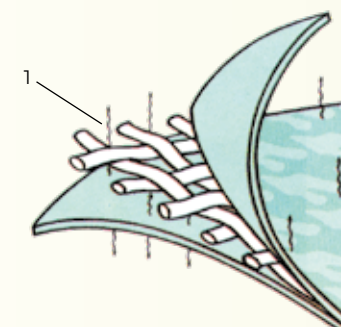

1. Sustainably integrated: active antibacterial components

Available in 5 colours



light blue white yellow beige mint



Cleanique®

CQ: ropimex®- seamless Hygienic-curtain
made of thermoplastic polyether polyurethane
for sanitary blocks and wet room

- Width: 190 cm
- Weight: 270 g/m²



- Care instructions



Available in 2 colours



velvet-white

sky blue

Standard dimensions

	Width	Height in cm
CQ108	100	80
CQ128	120	80
CQ121	140	120
CQ122	140	140
CQ123	140	170
CQ151	170	120
CQ152	170	140
CQ153	170	140
CQ211	220	120
CQ212	220	140
CQ213	220	170
CQ1220	120	190
CQ1420	140	190
CQ1620	160	190
CQ1820	180	190
CQ2020	200	190
CQ2220	220	190
CQ2420	240	190
CQ2720	270	190
CQ3020	300	190
CQ3320	330	190


Special dimensions on request.

ZEROTEX

ZEROTEX

ropimex® Zerøtex - non-woven
made of 100 % polypropylene (PP).
Saves time and money

- Width: ca. 210 cm
- Weight: 70 g/m²

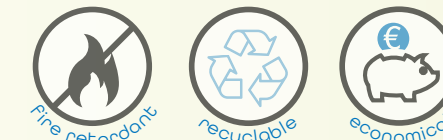

Suitable for
ropimex RTS telescopes and RTI as
well as the RVS rail system



Standard dimensions

	Width x Height in cm
	suitable for
Z1516	RTS 1.2 150 x 160
Z1517	RTI 1.3 150 x 175
Z1816	RTS 1.5 180 x 160
Z1817	RTI 1.6 180 x 175
Z2116	RTS 2.1 210 x 160
Z2117	RTI 2.2 210 x 175

Minimum order quantity 10 pieces /
per colour and size

Available in 5 colours



blue


yellow

light grey

apple green

creme

Folding Curtain ropicare

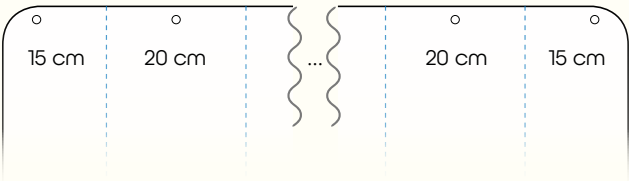

Slat width: 20 cm,
Start and end slat: 15 cm



Only in connection with ropimex® Vario RVS rail system, color coated white (see page 18 - 21)



ROC : ropimex®- Folding curtain made of PVC/polyester fabric. Easy care - no washing required



Scan

- Width: 180 cm
- Weight: 570 g/m²



Standard dimensions


	Width	Height in cm
ROC 110	110	175
ROC 130	130	175
ROC 150	150	175
ROC 170	170	175
ROC 190	190	175
ROC 210	210	175
ROC 230	230	175
ROC 250	250	175
ROC 270	270	175
ROC 290	290	175
ROC 310	310	175

Special dimensions on request.



Optionally available with a transparent area covering the entire width of the curtain.

Available in 8 colours

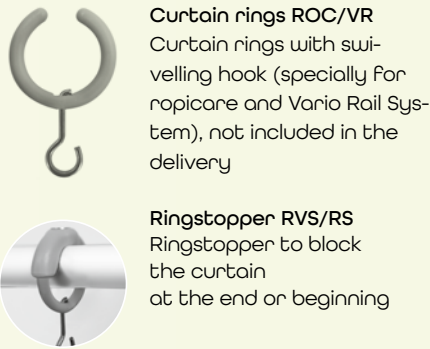

Curtain rings ROC/VR
Curtain rings with swivelling hook (specially for ropicare and Vario Rail System), not included in the delivery

Ringstopper RVS/RS
Ringstopper to block the curtain at the end or beginning

Folding screens RFW - wall-mounted



Scan

Infill panels:
non-see-through but
translucent

Non-see-through 360°
hinges

10 YEARS
warranty
on hinge
system

The robust and impact-resistant Folding screen is positioned in seconds and parked back against the wall after use.

Non-see-through but translucent leaf infill material with high-quality aluminium profile creates a pleasant, light atmosphere.

Modern leaf designs, available in different colours, give the folding screen an individual touch.



CE

Folding wall RFW wall-mounted



NEW!


Now with exclusive illustrations by Lisa Tegtmeier
The renowned Hamburg-based illustrator is known for her work for international companies such as Adobe, Aveda, Avon, Google, Klarna, Levi's, Nike ... and many more.



Folding screen for wall-mounting RFW CE

Leaf width	30 cm
Heights	145 / 165 / 185 cm
Weight	3,00 / 3,50 / 4,00 kg / leaf
Frame	Frame and hinged profiles: extruded aluminium (Al Mg Si 05), anodised surface, connecting screws: V2A stainless steel
Plastic parts	PA 6 GF 15, polyamide with 15% glass fibres
Infill panels	Polycarbonate, 2 mm, white, translucent, Flame resistant, construction material class B1, DIN 4102-01, extremely impact resistant
PFlege	commercially available cleaning product and disinfectant

Screen	open width	
3-leaf	90 cm	
4-leaf	120 cm	
5-leaf	150 cm	
6-leaf	180 cm	
7-leaf	210 cm	
8-leaf	240 cm	
9-leaf	270 cm	
10-leaf	300 cm	

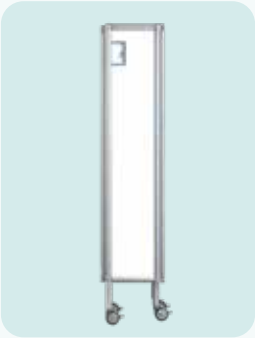

Non-see-through 360° hinges.
The RTC unit (ropimex®-Turn&Connect) enables 360° leaf rotation.

29 colours
23 designs
see pages 34-35



Folding screen RFW Accessories



Spacer set For wall fixing RFW-WA	Spacer For wall fixing, standard spacing 100 mm, special length to customer-specified dimensions possible up to 400 mm
---	--



Clamp set RFW-KLS	Clamp set for fixing to double-track standard rail
----------------------	--

Folding screen RFW column-mounted



Stainless steel column CE
For folding screen

Ideal wherever wall-mounting is not possible.
Stainless steel column with dimensions 5 x 5 cm.
Available in 3 heights.

RFW-ES4	Height 145 cm, weight 5,00 kg
RFW-ES6	Height 165 cm, weight 5,50 kg
RFW-ES8	Height 185 cm, weight 6,20 kg
Folding screen	Description see page 30



Stainless steel column with possibility to fix up to 4 Folding walls

Folding screens RFW - with trolley



Scan

Infill panels:
non-see-through but
translucent

Non-see-through 360°
hinges

10 YEARS
warranty
on hinge
system



Folding screen RFW with trolley: Position, Fold open, Fold together after use and roll away – care of the easy-clean surfaces made of aluminium and scratch-resistant polycarbonate is just as easy as the handling.



Trolley for Folding screen CE

The mobile version with max. 10 leaves enables mobile use. Please note: mobile folding screens are available in heights 145 and 165 cm only

RFW-FG4	Height 145 cm, weight 14 kg
RFW-FG6	Height 165 cm, weight 15 kg
Trolley base	52 x 59 cm
Folding screen	For description and screen length see page 30

Note:
Mounted on a
trolley, the folding
screen is only mobile.
It is therefore not
possible to convert
it to a wall-mounted
screen or to fix it to
an equipment rail or
column.







29 colours
23 designs
see pages 34-35



For easy handling - ingenious details such as lockable castors and lifting mechanism



Folding wall RFW 29Colours

				
01 Golden yellow	02 Bath tuquoise	03 Clove red	04 Lightblue	05 Copper brown
				
06 Vespa yellow	07 Opal turquoise	08 Red grey	09 Paris blue	10 Melon red
				
11 Apricot creme	12 Clown green	13 Tulle violet	14 Madonna blue	15 Apricot orange
				
16 Milk star white	17 Pale mountain lake turquoise	18 Rosewood apricot	19 Blue Grey	20 Mohair mauve
				
21 Nile green	22 Marker green	23 Sahara light red	24 Bermuda blue	25 Soft greige
				Standard: white translucent (without extra charge)
26 Foggy grey	27 Pale pistachio	28 Pepper white	29 Vintage white	30 Standard white translucent

Customize your folding wall online:




Scan

<https://ropimex.com/en/configurator/>



Colours based on RAL D2. Printed colours may deviate for technical reasons.

Folding wall RFW 23designs

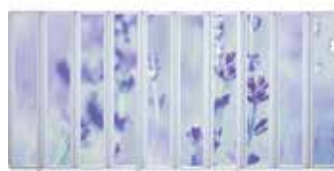structure			
	Brightwood	Darkwood	Abstract
Floral			
	Ginkgo	Yellowpalm	Dandelion
			
	Wildplants	Pivoine	Papillion
			
	Meadow	Lavender	Leaves
landscape			
	Relaxbeach	Heron	Town
			
	Flamingo	Bananaplants	Dreambeach
kids			
	Dino	Space	Submarine
Pop Art			
	Care	Clinics	

NEW!
Exclusive illustrations
by Lisa Tegtmeier

Butterfly Folding wall Model RBF

Infill panels:
non-see-through but
translucent

Non-see-through 360°
hinges

10
YEARS
warranty
on hinge
system



NEW!

Now with exclusive
illustrations by Lisa
Tegtmeier

The renowned
Hamburg- based
illustrator is known
for her work for in-
ternational compa-
nies such as Adobe,
Aveda, Avon,
Google, Klarna,
Levi's, Nike ... and
many more.




Scan

available in 2 dimensions

RBF



RBF200

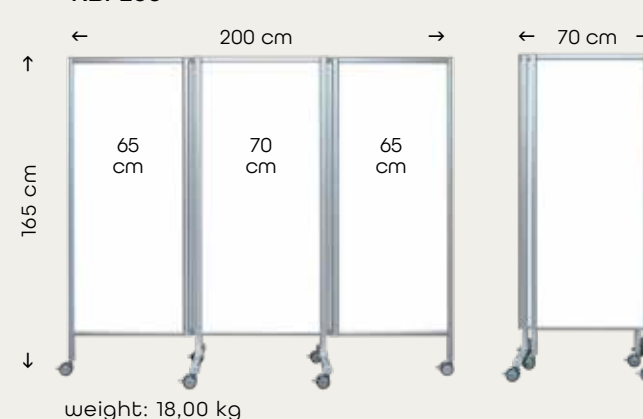

ropimex® Butterfly Folding wall RBF **CE**

3-leaf made in polycarbonate, with middle
leaf centred on double T-leg and 360°
rotatable leaves on either side, Fold-away.

Material	Frame and hinges: extruded alu- minium (Al Mg Si 05), anodised surface, connecting screws: V2A stainless steel, Infill panels Polycarbonate, 2 mm, Plastic parts PA 6 GF 15, polyamide with 15% glass fibres
Colour	Standard white, translucent. On option 29 colours, 21 designs



Butterfly Folding wall RBF / RBF 200 **29** colours

Colours and designs on option



Customize your Folding wall online:



Colours based on RAL D2. Printed colours may deviate for technical reasons.

Butterfly Folding wall RBF / RBF 200 **23** Designs



NEW!

Exclusive illustrations by Lisa Tegtmeier

Butterfly Paravent Model RBP

Non-see-through 360° hinges

10 YEARS warranty on hinge system



Lockable cas





Covering: staphCHEK



Scan

Available in 3 dimensions

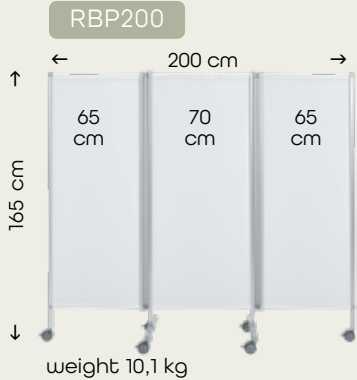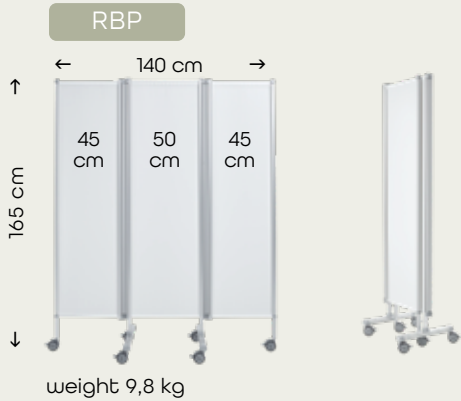


ropimex® butterfly paravent RBP CE

3-leaf, with middle leaf centred on double T-leg and 360° rotatable leaves on either side, Fold-away

Material	Frame: Aluminium design profiles, natural-anodised, made of polyamide, light grey On option: Plastic castors, Ø 65 mm lockable
Covering	Standard covering staph CHEK. staph CHEK Fabric is antibacterial, Flame-resistant, self-extinguishing, antistatic and durable. It does not have to be washed. Stains are easily wiped off with wipes. Available with ropitex-covering for an additional charge.
Colours	white, light blue, yellow, beige, mint

Recommended curtains see from page 22

VHC staph CHEK

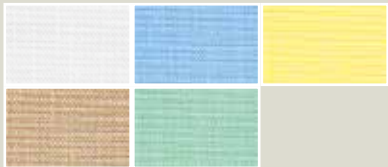


Lightweight paravent Model RLP

Non-see-through 360° hinges

10 YEARS warranty on hinge system



Paravent leaves must always be positioned at an angle to each other, unless the ropimex® base (model RLP/C) is used See page 45.

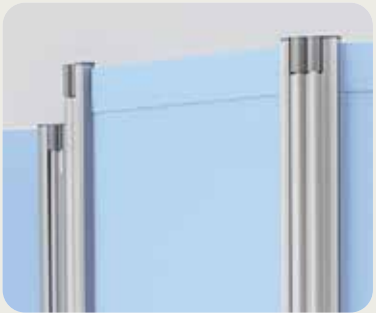

ropimex®-Lightweight paravents are convincingly easy to handle. The 360° non-see-through hinges each leaf enable each leaf to be moved into the required position. Covering without drapery.



ropimex® lightweight paravent RLP

Material	Frame: Aluminium design profiles, natural-anodised, Foot and cover caps made of polyamide, light grey. On option: Plastic castors, black	
Covering	Standard covering staph CHEK. staph CHEK Fabric is antibacterial, flame-resistant, self-extinguishing, antistatic and durable. It does not have to be washed. Stains are easily wiped off with wipes. Available with ropitex-covering for an additional charge.	Recommended Curtains see page 23
Colours	white, light blue, yellow, beige, mint	




Base to stabilise the lightweight paravent Model RLP/C



Base RLP/C



holding up to 4 paravents

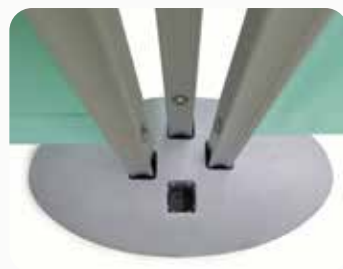

3 paravents in the base





4 paravents in the base



Paravents in a row as a routing system



ropimex® base for lightweight paravents

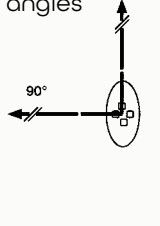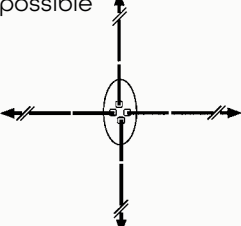Model RLP/C	Base to stabilise the paravent when opened out straight (use of the open width) and/or for holding several paravents (max. 4 per base), as a partition or for setting up a routing system		
Material	Glass-fibre reinforced polyamide (GRP)		
Dimensions/Weight	8 x 15 x 30 cm (H x W x D), weight 1,3 kg		
Variations	1 Paravent in a line  1 Paravent at right-angles  up to 4 paravents possible 		

Lightweight paravent for wall-mounting Model RLP/W



ropimex® Lightweight paravent RLP with wall connector

Model	Leaf	Weight	Width	Height
RLP/W2	2	3,8 kg	100 cm	165 cm
RLP/W3	3	5,7 kg	150 cm	
RLP/W4	4	7,6 kg	200 cm	
RLP/W5	5	9,5 kg	250 cm	

Accessories

RLP-WA	Spacer for wall fixing, standard spacing 100 mm, special length to customer-specified dimensions possible up to 400 mm	
RLP-KLS	Clamp set for fixing to double-track standard rail	
RLP-ES	Upright made of stainless steel for RLP/W, for free-standing installation in the room, height 1,650 mm	

Shower splash guard
Telescopic screen
cubicles

Model **RDSE/RDLE/RDTE**

Model **RKLE**

Wall mounts
WH/F, WH 95 see p.15

Curtain rings
included in the scope of
delivery. Colour grey.

Telescopic screen cubicles

With this space-saving model it takes just a few seconds to create a screen for washing, changing either as a corner solution or as a U-shaped solution.

Model	Corner solution	
		recommended wall mount
RKLE 1212	Telescopic screen cubicle, max. 120/115 cm, min. 75/70 cm	WH 95 or WH/F
RKLE 1515	Telescopic screen cubicle, max. 150/145 cm, min. 90/85 cm	
	U-shaped solution	
		recommended wall mount
RKAE 1.2	Telescopic screen in combination with RKLE 1212 as a U-shaped solution 120x115x120 cm	2 x WH/F
or RKAE 1.5	lescopic screen in combination with RKLE 1515 as a U-shaped solution150x145x150 cm	



Shower splash guard

Model	Telescopic extension		
RDSE	Splash guard arm swing-out and Fold-down, dimension 100 cm		recommended wall mounts
RDTE	Splash guard arm swing-out, dimension 95 cm		
RDLE	Splash guard arm, L-shaped, can be used on the right and left-hand side, side leg 25 cm, swing-out and Fold-down, dimension 95 cm		WH 95 or WH/F




Recommended curtains for wet room areas: CQ, VH/C see p. 23 -24

CQ Cleanique VHC staph CHEK



Recommended curtains for changing rooms: RTX see p. 22

RTX ropitex

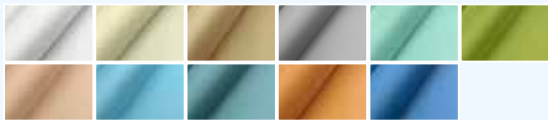

Curtains washable

Surfaces disinfectable

with steam disinfectable

with vaporison disinfectable

sustainable repairable

ropimex® convincing product quality for maximum customer satisfaction

ropimex® HyGentle: Product design and material combination for maximum hygiene

Product design: Hygiene starts with the design. Corners, angles or geometries that are ill-suited for cleaning and disinfection are avoided to the greatest possible extent.

Material composition: Cleanability and disinfectability are prioritized in the selection of the materials required for manufacturing the product. This is aimed at achieving an optimum where soiling and microbiological pollution are concerned. Porous, absorbent surfaces must be avoided.

Cleaning and disinfection: We test and document the materials of our products in terms of their natural microbial contamination, just as we document their cleanability and disinfectability with commercially available detergents and disinfectants.

ropimex® business divisions at the Neunkirchen site




ropimex® Screen Systems
privacy protection systems and flexible partitions
for healthcare facilities



ropimex® Sound Protection
Privacy and sound insulation systems
for the office area



ropimex® Chemicals
Disinfectants, refinements, Mosquito and
tick protection

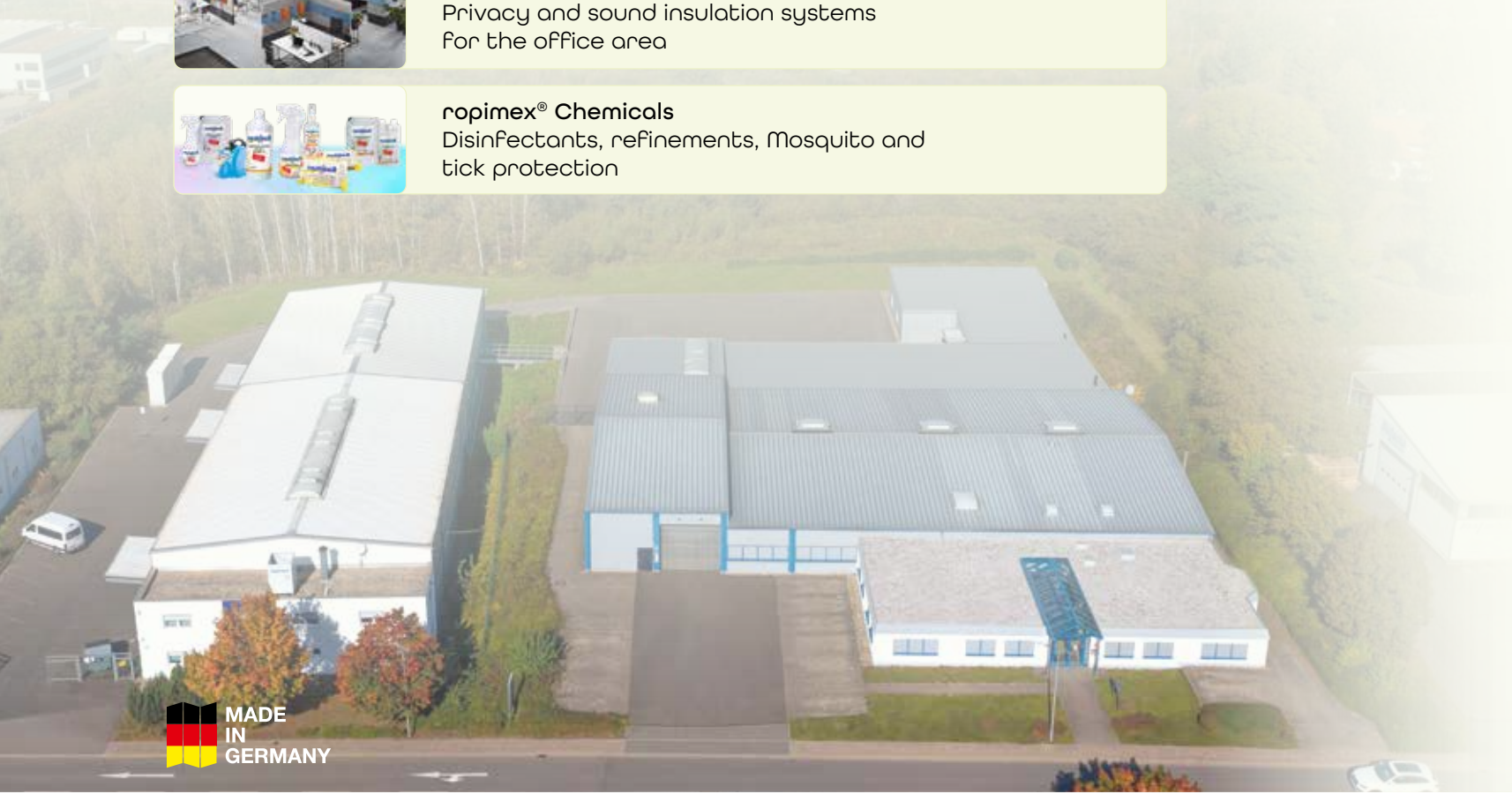

Satisfied customers around the globe



Since 1976, ropimex® has been busy finding solutions to protect the privacy of the sick and those in need of care.

The telescopic screens of the early years are still in use today. From them a system developed whose diversity, use and combination options provide a suitable solution for each and every customer.

All these properties, top advice from authorised stockists and the fastest possible delivery have secured ropimex® satisfied users in over 60 countries, all over the world.



ropimex® Our commitment



Scan



Ecological / sustainable

- ropimex® products are 98 % sourced in Germany and 100 % manufactured at the plant in Neunkirchen/Saar.
- consistent recycling of our waste such as aluminium, fabric scraps, cardboard and plastic
- we are greatly reducing our paper consumption. Our goal: paperless office
- gradual conversion of electricity-intensive consumers
- our products are phthalate-free



Tree sponsor - ropimex® supports reforestation project from PLANT-MY-TREE®



Social behaviour

- Well-being is a priority. We are happy to invest in the environment of the company
 - whether it's spacious, open-plan offices with modern equipment,
 - friendly communal areas,
 - flexible working hours or
 - homeoffice
- thanks to these practical measures, the sickness rate is low, the fluctuation rate is virtually zero.



Commitment

We also want those to benefit who are unable to participate in the performance process, or only to a limited extent. For years, we have been involved as a partner in the Lebenshilfe organisation and other measures, particularly those that support young people.

- Partner-Company 2024 - Ropimex R. OPEL GmbH is a Partner-Company of DOCTORS WITHOUT BORDERS / MEDECINS SANS FRONTIERES, we support their work with our donation.

DOCTORS WITHOUT BORDERS provides emergency medical relief in areas fraught by crisis, war, or natural catastrophes. The international organisation helps people quickly, effectively and unbureaucratically, irrespective of their race, religion or political orientation. For further information: www.msf.org



Medical device

ropimex® products meet the requirements of the regulatory specifications for medical devices in accordance with the annex to Declaration of Conformity EU 2017/745, ropimex® visual protection systems.

The systems are designed so that a patient zone / intimate zone can be created for virtually any existing condition - stationary or mobile.

ropimex® screen systems are customisable accessories for medical devices such as patient beds, treatment tables and treatment chairs in medical facilities.



Certified according to:


DIN EN ISO 9001:2015 and DIN EN ISO 14001:2015

The certification body of TÜV SÜD Management Service GmbH certifies that the company ROPIMEX R. OPEL GmbH, Bildstocker Str. 12-14, 66538 Neunkirchen, Germany has **introduced and applies a quality management system** for the scope of manufacture and sale of privacy and soundproofing elements. An audit was performed, Order No. 707145152. Proof has been furnished that the requirements according to DIN EN ISO 9001:2015 are fulfilled.

Furthermore, the company has **introduced and applied an environmental management system** for the production and sale of privacy screens and soundproofing elements as well as biocidal products, disinfectants and coating systems. An audit has provided proof that the requirements of DIN EN ISO 14001:2015 have been met.

For your safety

ropimex® products meet the highest safety standards.

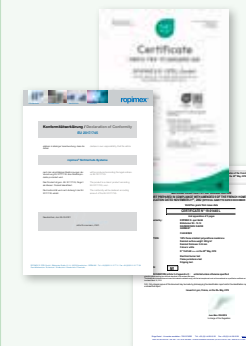

Screen- Systems

Healthcare | 2025



Telescopic screens | Telescopic arms | Mobile screens | Vario rail system
Curtains | Folding screens | Paravents | Shower splash guard



Scan

ROPIMEX R. OPEL GmbH
Bildstocker Str. 12-14 · D-66538 Neunkirchen
Tel. +49 6821 9 12 77-0
info@ropimex.com · www.ropimex.com

ropimex[®]
since 1976