

ScanAfia XL nursing bed

A Quality New Arrival for Home and Home-like Care



Main features:

1. Mattress base width 120 cm for increased comfort and safety
2. Very high Safe Working Load (SWL) of 285 kg enables care of larger patients
3. Impressive height adjustment range improves patient independence and carer comfort

The special sized ScanAfia XL nursing bed has been designed with larger and heavier patients in mind. The bed is particularly suitable for nursing homes, where the average size of patients has increased. The special sized bed is 120 cm wide, with a maximum patient weight of 250 kg.

The ScanAfia nursing bed can be combined with our range of tables and chairs, creating a harmonious décor without compromising functionality. The solid white oak bed ends give the bed a longer life cycle than standard beds. Our high quality epoxy coating also contributes to the long lifespan of Lojer beds. The bed features a safe 24V low voltage electrical system, and its electrical motors and control units are splash resistant. The design of the bed is based on the 60601-1-52 nursing bed standard. The large twin castor make the bed easy to manoeuvre yet extremely sturdy when locked in position.

LOJER® *For easy care*


Technical data and features	
Overall width with rails	130 cm
Overall length	217 cm
Height adjustment range	34-82 cm (150 mm twin castors) 32-80 cm (optinal 125 mm twin castors)
Recommended mattress size	205 x 118 x 13 cm
Total weight (inc. rails)	145 kg
Castors	Ø150 mm central locking twin castors and direction wheel
Safe Working Load (SWL)	285 kg
Back section adjustment	70°
Thigh section adjustment	30°
Calf section adjustment	11°
Other standard features	Solid wood bed-ends of white-oak double-locking sides with locking button (pair) hand control manual lowering of back section anti-Trendelenburg integrated mounting for lifting pole and IV pole
Electrical system	24V
CE-marked	Yes
Warranty	2 years
Country of Manufacture	Finland

Factory installed options:

- Battery for back-up
- Ø125 mm twin castors
- Set of bumper rollers

Other accessories:

- Lifting pole with triangle handle
- Lift support
- IV pole with one hand quick adjustment
- Reading lamp
- Extra height side rail
- Padded coating on top of the side rail, one unit
- Urine bag holder
- Nurse locking panel (allows functions to be locked to prevent patient use)
- Flexible hand controller holder

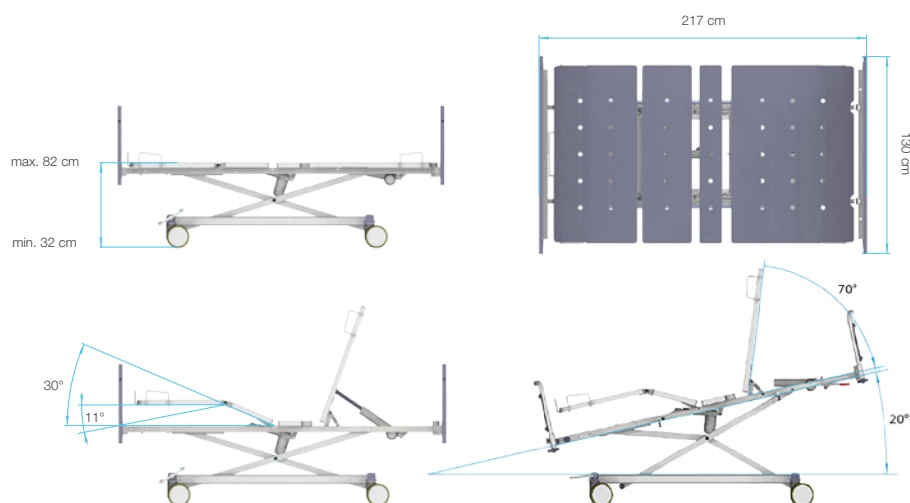

LOJER® For easy care